


ELKEINE VR PH-T5511/5521



3 Traction modes

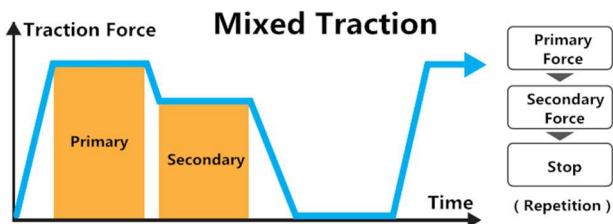

The traction mode consists of 3 tractions which are consecutive traction, intermittent traction and mixed traction. The setting content is different which depends on the mode.

Easy to use control panel



Simple designed control panel makes the operation much easier. 8 kinds of ending sounds help to identify easily, even when multiple traction force devices are being used.

Cold feeling relieving function



Bed built-in type "Cold feeling relieving function" is standardly equipped. Because of the built-in type, there is no need to attach it every time and it will not make the lumbar feel uneven. In addition, "Cold feeling relieving function" which is used for warming the neck is also standardly equipped.

One-touch underarm brace



It is designed to pursue both the functionality and safety. The patients can grasp the brace during the treatment which makes it more steadily.

Cervical vertebrae orthosis & lumbar brace



EXGEL® which has characteristics of body pressure dispersion, shape retention and shock absorption is used to the lower jaw part. It makes the jaw much fitter to the cervical vertebrae orthosis and prevents too much pressure on the lower jaw.



Anti-skid tape is attached inside of the lumbar brace which is used to fix the waist, which makes the waist fixed firmly and difficult to slip off.

Safety function



Control circuit provides multiple safety protection from all angles. When the abnormality is detected, it works when necessary and automatically controls the traction force.



Besides the control circuit, the emergency stop switch which is located on the top cover can also be manually pushed by the staff to stop the actions of the device.

## Optional color variation for different rehabilitation rooms

For accessories of “flatbed”, “one-touch underarm brace”, “traction triangular pillow”, and “traction pillow”, besides the standard color “purple”, 3 option color can also be chosen. (The color of “cervical vertebrae orthosis”, “regular side brace”, “lumbar brace”, “head protection belt”, and “cervical traction chair” is gray.)



Light green



Orange



Pink



The color of the upper parts is gray.

## Optional layout according to the treatment content and rehabilitation rooms

ELKEINE VR is designed to meet many different layout variations. The most suitable type can be chosen according to the treatment content.

### ELKEINE VR PH-T5511 (1 CH, one-person use)



Lumbar / Cervical (A Type)



Lumbar / Cervical (AB Type)



Cervicalx1 (KS Type)



Lumbarx1 (C Type)

### ELKEINE VR PH-T5521 (2 CH, two-person use)



Lumbar & Cervical (A Type)



Lumbar & Cervical (AB Type)



Cervicalx2 (K Type)



Lumbarx2 (WHA Type)

## Specification

	ELKEINE VR 5511	ELKEINE VR 5521		
Commodity code	PH-T5511 (100V) PH-T5511EK (220V)	PH-T5521 (100V) PH-T5521EK (220V)	<b>Traction force</b> 9.8-970N=1-99kgf (9.8N=1kgf) ±20%	
			<b>Number of channels</b> 1 CH (One-person use)	2 CH (Two-person use)
Main power	AC100V or 220V 50/60Hz		<b>Treatment mode</b> Intermittent, mixed, consecutive	
Power consumption	105VA (100V) 100VA (220V)	190VA (100V) 180VA (220V)	<b>Treatment time</b> 1-30 minutes(1 minute unit) or consecutive ±10%	
	Main unit: 574(W)×304(D)×860(H)mm Flatbed: 580(W)×2,015(D)×555(H)mm		<b>Electrical class, type</b> Class I, Type B	
Weight	Main unit: 55kg Flatbed: 41kg	Main unit: 65kg Flatbed: 41kg	<b>Safety function</b> Abnormal operation detection circuit, motor detection control circuit, patient stop switch, emergency stop switch	



**IMPORTANT** \*Please read the User Manual carefully before using our products.  
**NOTICE** \*The specifications and design may change without notice.

Made in Japan

WKRF-PHT5500-IS

IS-1711-00 I

**NIHONMEDIX**

<http://www.nihonmedix.co.jp/en/>